

Project: Sapien 190
Date: 16.05.2010
Version: 0.05

Product architecture document for Sapien 190

Embedded Real-Time Systems (TI-IRTS) Spring 2010

Peter Høgh Mikkelsen (20087291)
Anders Block Arnfast (20085515)
Kim Bjerge (20097553)

PAsien

Document history

Date	Version	Description	Author
24.04.2010	0.00	Initial Version	KBE
10.05.2010	0.01	Updated with Use Case View and initial UML diagrams for first delivery UC#1	KBE
15.05.2010	0.02	Updated diagrams with UC#1 and added Ipump protocol in □omm.. Package	KBE
15.05.2010	0.03	Added command pattern for setting parameters	KBE
16.05.2010	0.04	Started on chapter 13 and 14	KBE
16.05.2010	0.05	Added details for chapter 13 and 14. Still missing error messages.	KBE

Template User guide:

- To be deleted in the final document

This Template is a highly modified version of a RUP template (Rational Unified Process) for a software architecture document. The modifications includes elements from the Danish Structured Development Handbooks (SPU) SW documentation guide. The template is based on the famous “4+1” View Model described by Philippe Kruchten, “The 4+1 View Model of Architecture”, IEEE Software Nov.1995.

The Template is written as a **Microsoft Word .dot** file. A double click with the mouse on the xx.dot file in the file manager will open a new word document with the default name: Document 1, including the information in the template.

Start with modifying the document title and subject. Document title and subject (e.g. System/product <name>) is modified in **File/Properties – Summary**. Title and subject will automatically be inserted on the front page and in the page header. Save the document with a new name.

If the **.dot** file is installed in the directory where word looks for template files, it will be possible to activate a new document based on this template by selecting **File/New** from the Word program.

For helping the user with writing the document, help text is included for each section in the document. This help text is written as “Hidden Text” in Word.

To show (enable) the help text on the screen:

Menu choice: **Tools/Options – View**

and check the check box: Hidden text (alternatively remove the mark – for removing the help text).

To print the template file with help text:

Menu choice: **Tools/Options – Print**

and check the check box Hidden text.

To modify the template:

Start Word and open the template by selecting: **File/Open**

Table of Contents

1. INTRODUCTION.....	4
1.1 Purpose and Scope	4
1.2 References	4
1.3 Definitions and acronyms	4
1.4 Document structure and reading guide	5
1.5 Document role in an iterative development process	5
2. SYSTEM OVEVIEW	5
2.1 System context	5
2.2 System introduction	5
3. SYSTEMET INTERFACES	6
3.1 Interface to human actors	6
3.2 Interface to external system actors	6
3.3 Interface to hardware actors	6
3.3.1 Pulse interface	6
3.3.2 Infusion pump interface	6
3.4 Interface to external software actors	6
4. USE CASE VIEW	6
4.1 Overview of architecture significant Use cases	7
4.2 Use case #1: Execute and Control Simulation scenarios	7
4.2.1 Use case goal	7
4.2.2 Use case scenarios	7
4.3 Use case # 2: Select and Initiate Scenario scenarios.....	8
5. LOGICAL VIEW.....	9
5.1 Overview	9
5.2 Architecturally significant design packages.....	10
5.2.1 Continuous Package.....	10
5.2.2 Discrete	14
5.2.3 Com.....	15
5.2.4 AbstractHW	16
5.2.5 AbstractOS	16
5.2.6 Application Helper Classes.....	17
5.3 Use case realizations	18
5.3.1 Use case #1: Execute and Control Simulation scenarios	18
5.3.2 Use case 2. realization	24
6. PROCESS/TASK VIEW	25
6.1 Process/task overview	25
6.2 Process/task implementation	25
6.3 Process/task communication and synchronization.....	25
6.4 Process group 1.	25
6.4.1 Process communication in group 1	25
6.4.2 Process 1. description	25
6.4.3 Process 2. description	25
6.5 Process group 2.	26

7. DEPLOYMENT VIEW	26
7.1 System configurations overview	26
7.2 System configurations	26
7.2.1 Configuration 1.	26
7.2.2 Configuration 2.	26
7.3 Node descriptions	26
7.3.1 Node 1. description	26
7.3.2 Node 2. description	26
8. IMPLEMENTATION VIEW	27
8.1 Overview	27
8.2 Component descriptions	27
8.2.1 Component 1	27
8.2.2 Component 2	27
9. DATA VIEW	27
9.1 Data model	27
9.2 Implementation of persistence	27
10. GENEREL DESIGN DECISIONS	27
10.1 Architectural goals and constraints	27
10.2 Architectural patterns	27
10.3 General user interface design rules	28
10.4 Exception and error handling	28
10.5 Implementation languages and tools	28
10.6 Implementation libraries	28
11. SIZE AND PERFORMANCE	28
12. QUALITY	28
13. COMPILATION AND LINKING	28
13.1 Rhapsody modeling and testing	28
13.2 Linux host Compilation-software	30
13.3 Linux Cross Compilation and linking process	30
13.3.1 Qt Cross Compilation with qt-everywhere	31
13.3.2 WFDB Installation and Cross Compilation	32
14. INSTALLATION AND EXECUTING	34
14.1 Installation	34
14.2 Executing-hardware	34
14.3 Executing-software	34
14.4 Execution-control (start, stop and restart)	36
14.5 Error messages	36
15. APPENDICES	36

1. INTRODUCTION

1.1 Purpose and Scope

The Sapient 190 is a patient simulator, simulating human physiological behaviour according to different patient scenarios. Scenarios are written in an open format and can be downloaded to the Sapient 190 that contains patient records from the PhysioBank database.

In conjunction with the Sapient 110 human doll, the Sapient system provides a complete simulated human interface, with EKG, ECG and respiratory measuring spots and medicine injection spots.

1.2 References

- [1] Erich Gamma et al., Design Patterns: Elements of Reusable Object-Oriented Software, Addison Wesley (GoF)
- [2] Bruce Powell Douglass, Real-Time Design Patterns: Robust Scalable Architecture for Real-Time Systems
- [3] PhysioNet and PhysioBank the research resource for complex physiologic signals
<http://www.physionet.org/>
- [4] Project specification for PSIMU: Patient Simulator System
<http://kurser.iha.dk/eit/tiirts/Projekter/PSIMU-project.doc>
- [5] Project specification for LMON: Local Monitor System
<http://kurser.iha.dk/eit/tiirts/Projekter/LMON-project.doc>
- [6] Project specification for IPUMP: Infusion Pump System
<http://kurser.iha.dk/eit/tiirts/Projekter/IPUMP-project.doc>
- [7] Project Interface Specification
<http://kurser.iha.dk/eit/tiirts/Projekter/ProjectInterfaces.doc>
- [8] Requirement specification for Sapient 190
- [9] Getting started with Qt
<http://devkit8000.wikispaces.com/Qt+and+Qt+Everywhere>

1.3 Definitions and acronyms

- **Model** – A model that represents one physical individual. The model may consist of sub-models for different body subsystems. The model uses an algorithm to compute the output signals based on the input signals or patient records.
- **Model Parameters** – Parameters used in the model.
- **Record** – A patient record file taken from the PhysioBank database
- **Scenario** – Model parameters and a collection of records taken from the PhysioBank database.
- **Scenario Configuration** – A set of files that represents model parameters and patient records.
- **Signals** – Signals in form of waveform files or input from external equipment
- **Simulation** – A continuous mode, where the physiological output signals are updated according to the model and the scenario applied to it.
- **EKG** – Electrocardiogram
- **EDR** – ECG-Derived Respiration

1.4 Document structure and reading guide

Chapter 2 and 3 gives an overview of the product and the interfaces to the patient simulator in terms of other devices and user operation.

The document describes the design using the “4+1” view. For each iteration as specified in ROPES [2] we have selected one or more use cases that is used in describing the Use Case, Logical, Process, Deployment and Implementations Views. These views are described in chapters 5 – 8.


1.5 Document role in an iterative development process

This document is updated in using the ROPES development process [2]. The document is updated for every design cycle of the ROPES spiral microcycle. Every design microcycle is using the 4+1 view

2. SYSTEM OVEVIEW

The Sapien 190 is a compact unit with a graphical user interface (GUI) and interfaces to a range of body monitoring and injection equipment. The Sapien 190 can be operated by the instructor and connected to a bedside monitor that used the output signals from the patient simulator. Optional a medicine infusion pump can be connected to the simulator.

2.1 System context

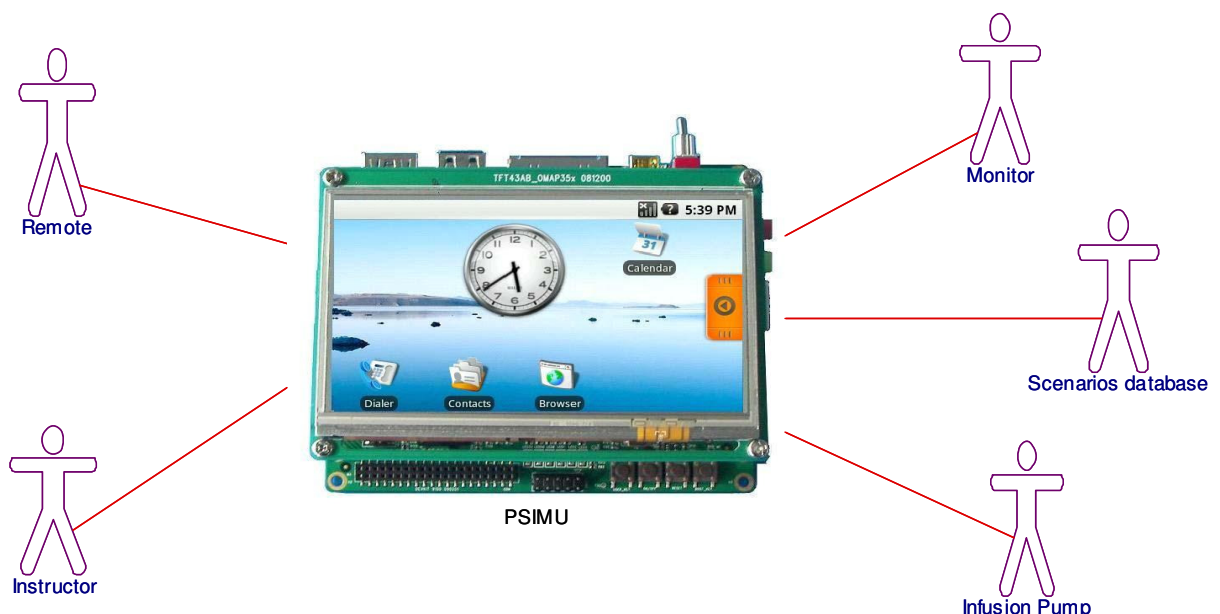

2.2 System introduction

The instructor will be able to control the simulation like forcing a heart attack and monitor the waveform of the signals that is send to the bedside monitor. Two analogue outputs signals can be connected to a bedside monitor to display ECG and EDR signals. It is possible remotely using an Ethernet connection to the patient simulator to update and delete files in the scenarios database on Sapien 190. The optional infusion pump injects medicine into the patient and the flow of medicine can be monitored on the Sapien 190 LCD display.

3. SYSTEMET INTERFACES

3.1 Interface to human actors

Instructor controls the patient simulation by inputs to the LCD touch screen....

3.2 Interface to external system actors

Remotely it is possible to update the scenarios database by use of a ftp connection with a Ethernet connected to the patient simulator.

3.3 Interface to hardware actors

Analogue outputs signals to the LMON has a resolutions of 12 bits sampled at 250 hz.
The analogue output has a voltage range of 0 - 4.096 volts.

3.3.1 Pulse interface

The pulse is send to a connected monitor using a RS232 connection (115200 baud, 8 databits, even parity, 1 stop bit).

Format of the pulse signal is a maximum 3 bytes ASCII value (beats/minute). It is transmitted from the PSIMU to the LMON with a rate of 1 second. The pulse value is terminated with a <CR>. <Pulse> <CR>

Example a terminal can be connected to display the pulse values updated each second:

103

107

3.3.2 Infusion pump interface

The infusion pump is connected using a RS232 connection with the protocol described below.

A fixed 64 bytes ASCII PDU is transmitted from IPUMP to PSIMU.

This PDU is transmitted every second when the pump is started.

PDU format:

4 bytes startframe (##?*)

46 bytes medicin name

12 bytes volume infused (since started)

2 bytes CRC checksum

3.4 Interface to external software actors

4. USE CASE VIEW

In the first iteration the UC #1 has been selected since it provides the main functionality and is the essential use case for the architecture of the patient simulator.

Execute and Control Simulation

When the patient simulation is running it will perform reading of the digitized physiologic signals and send these “real time” values as analogue signals to local connected bedside monitoring equipments. Up to 2 analogue channels with different signals is possible to be simulated simultaneously. The simulated signals can be ECG or EDR. The pulse will be signaled to the bedside monitoring equipment as a digital signal.

4.1 Overview of architecture significant Use cases

The UC #1 has been selected for the first iteration of the ROPES spiral microcycle [2] in making the architectural, mechanistic and detailed design. This use case is significant and provides the central functionality of the patient simulator. This use case provides the basis functionality that allows the monitor to be connected being able to display the ECG and EDR signals. It also reduces the risk for developing the patient monitor since it covers all the unknown technologies of the product like:

- Reading the patient record files on the target (Linux, WFDB and target)
- Generating the analogue output signals (Writing to drivers)
- Display of signal waveform using Qt on target (Working with Qt on target)

4.2 Use case #1: Execute and Control Simulation scenarios

4.2.1 Use case goal

To simulate signals from the patient based on the selected scenario.
The waveform of signals to monitor and the status of patient are displayed.
The instructor must be able to monitor the simulated patient.

4.2.2 Use case scenarios

Scenario 1 - normal:

1. Opens scenario file
2. Search for record file
3. Opens record file
4. Continues to reads samples from record file and performs:
 - a. Use patient model to generate ECG signal
 - b. Use patient model and ECG signal to calculate EDR signal
 - c. Use patient model and ECG signal to calculate pulse
 - d. Update output signals: Pulse, ECG and EDR

Scenario 2 – normal with alternative scenario record:

1. Opens scenario file
2. Search for record file
3. Opens record file
4. Continues to reads samples from record file and performs:
 - a. Use patient model to generate ECG signal
 - b. Use patient model and ECG signal to calculate EDR signal
 - c. Use patient model and ECG signal to calculate pulse

- d. Update output signals: Pulse, ECG and EDR
5. Open alternative record file after specified simulation time and continues at 2.

Scenario 3 – normal with IPUMP:

1. Opens scenario file
2. Search for record file
3. Opens record file
4. Continues to reads samples from record file and performs:
 - a. Reads medicine and volume from IPUMP
 - b. Use patient model to generate ECG signal
 - c. Use patient model and ECG signal to calculate EDR signal
 - d. Use patient model and ECG signal to calculate pulse
 - e. Update output signals: Pulse, ECG and EDR

Scenario 4 – error opening record file:

1. Opens scenario file
2. Failed to search and open record file
3. Error message on LCD display


4.3 Use case # 2: Select and Initiate Scenario scenarios

Initiate the simulation using a scenario configuration chosen by the instructor.

This use case is not part of the first iteration and will be updated in the next iteration.

5. LOGICAL VIEW

5.1 Overview

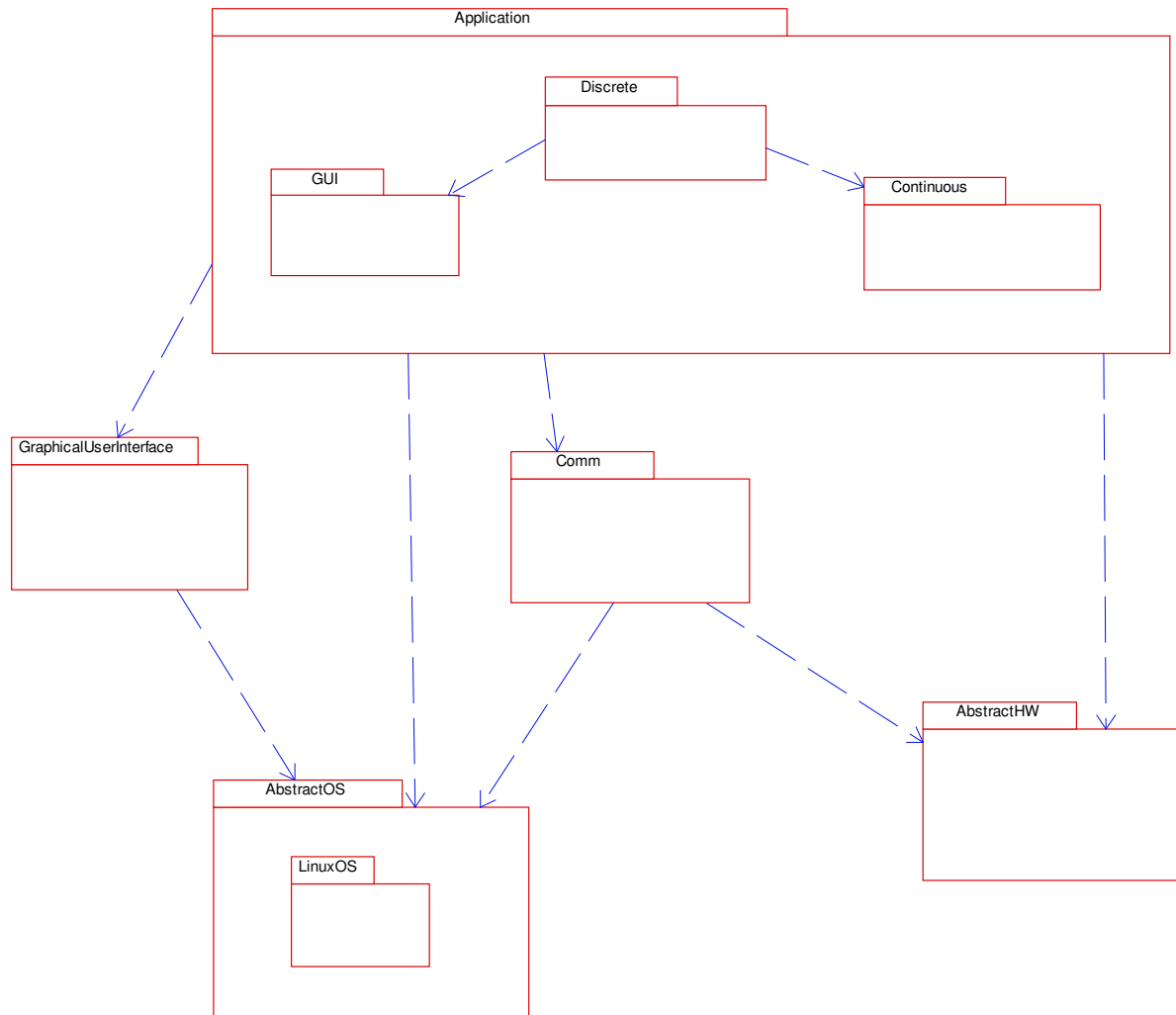

5.2 Architecturally significant design packages

5.2.1 Continuous Package

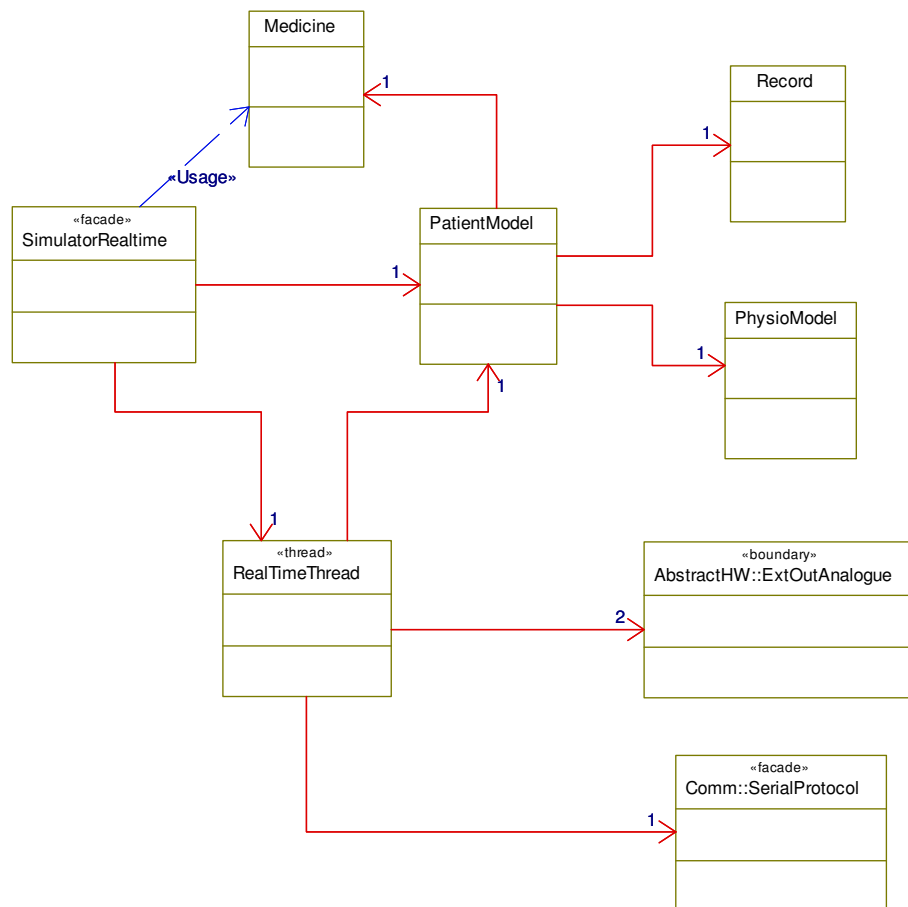


Figure 1 Major Classes in Continuous Package

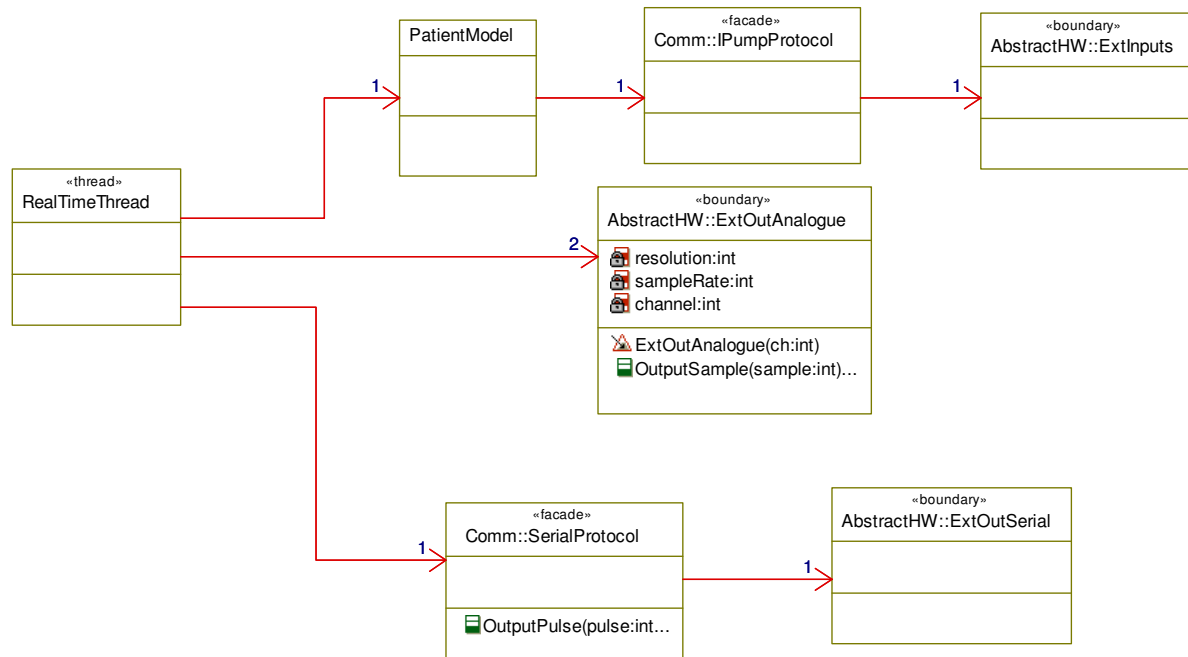


Figure 2 Boundary Classes for RealTimeThread

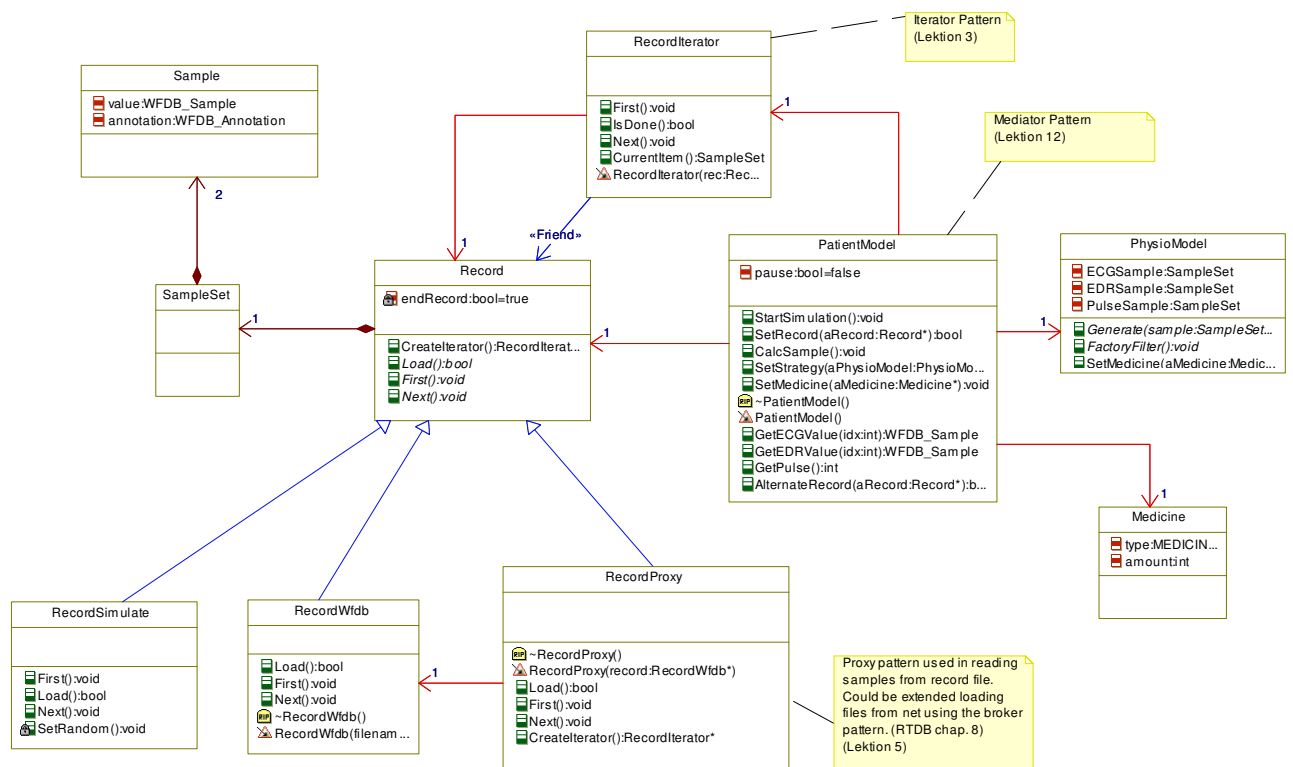


Figure 3 Proxy and Iterator Pattern used to access Records from the PatientModel

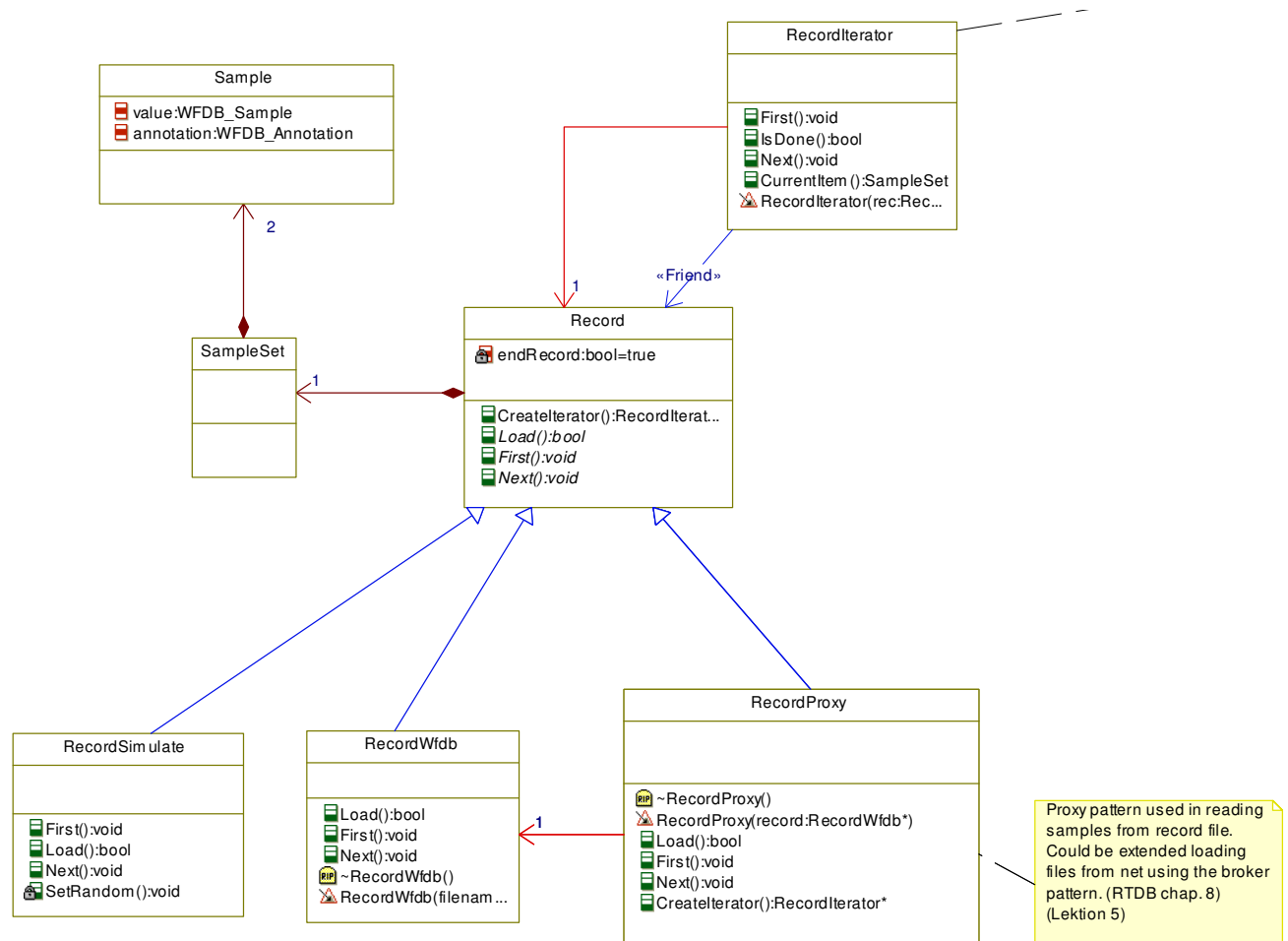


Figure 4 Record and Iterator

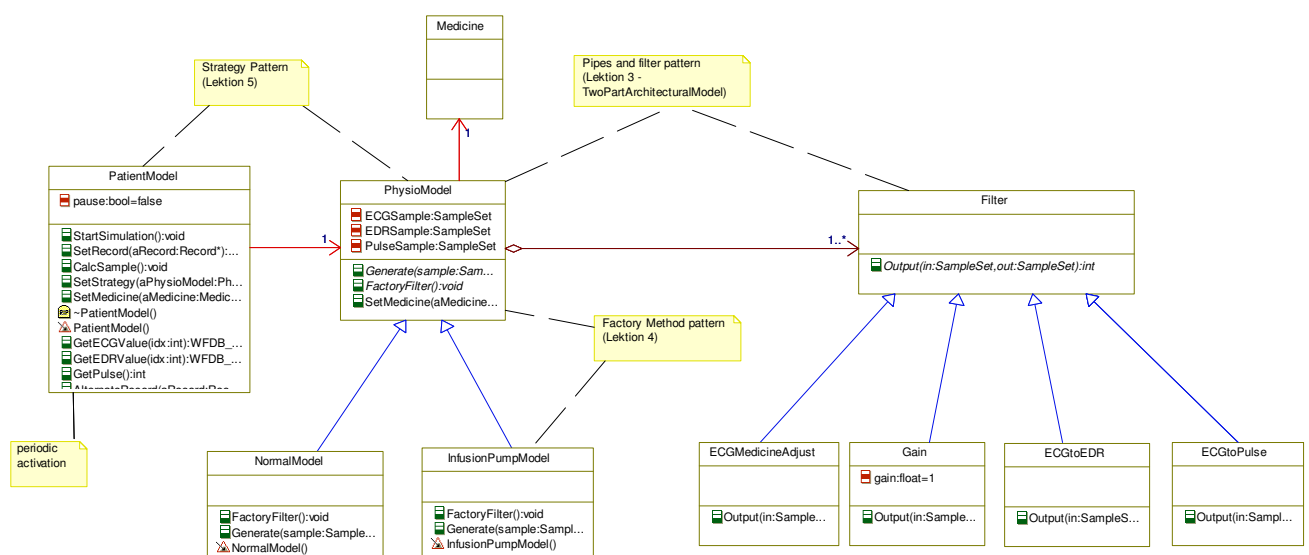


Figure 5 Strategy, Filter and Pipes Pattern used for PhysioModel

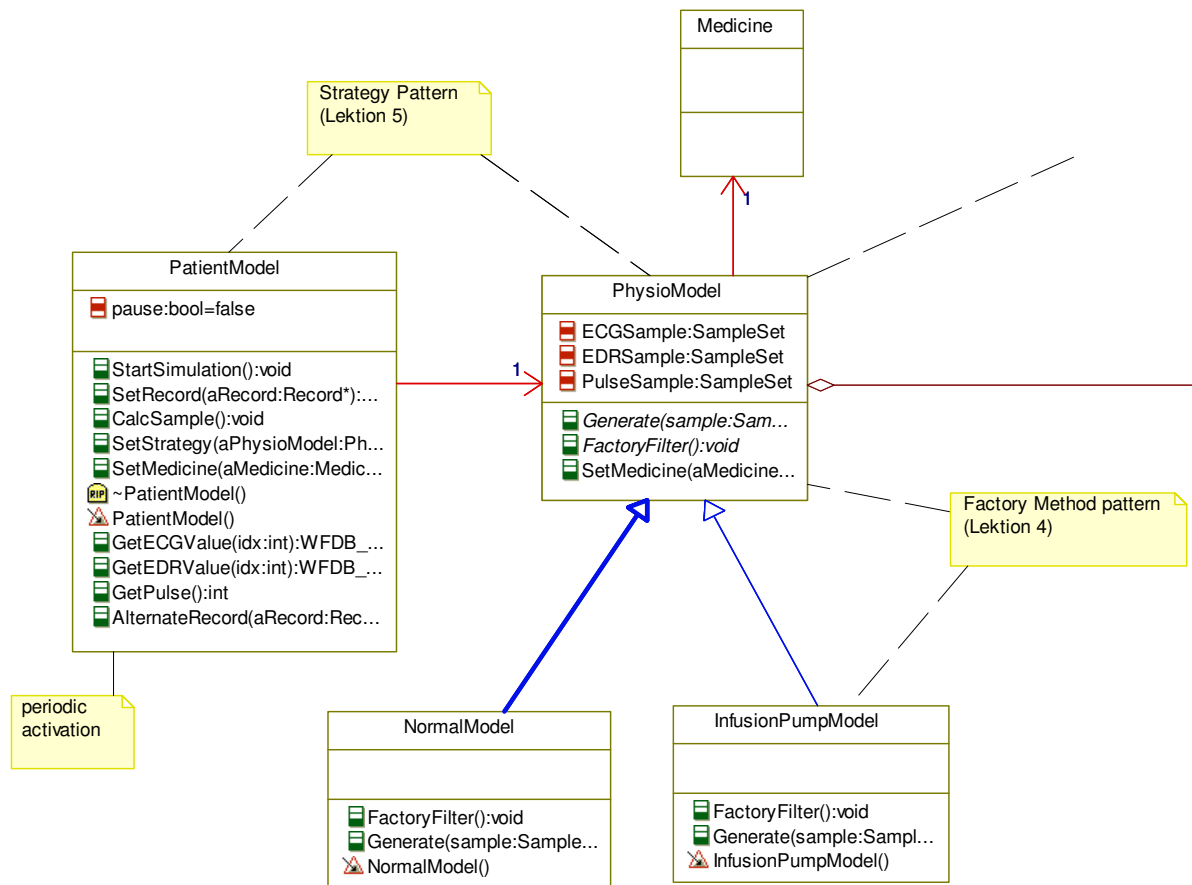


Figure 6 Strategy for PatientModel

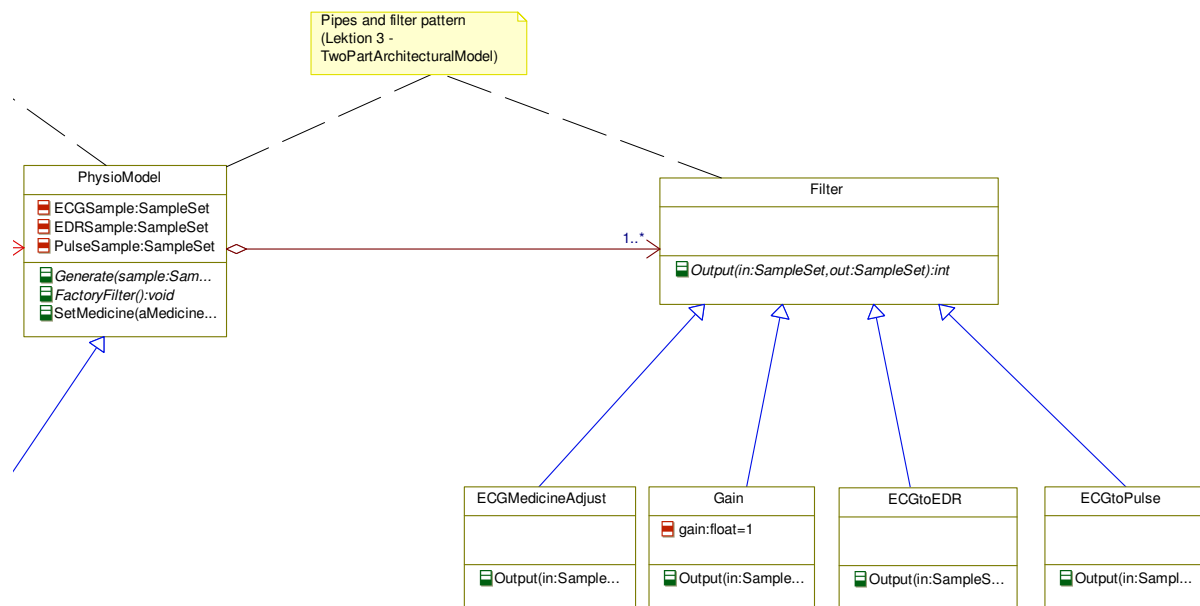


Figure 7 Filter Pattern for PhysioModel

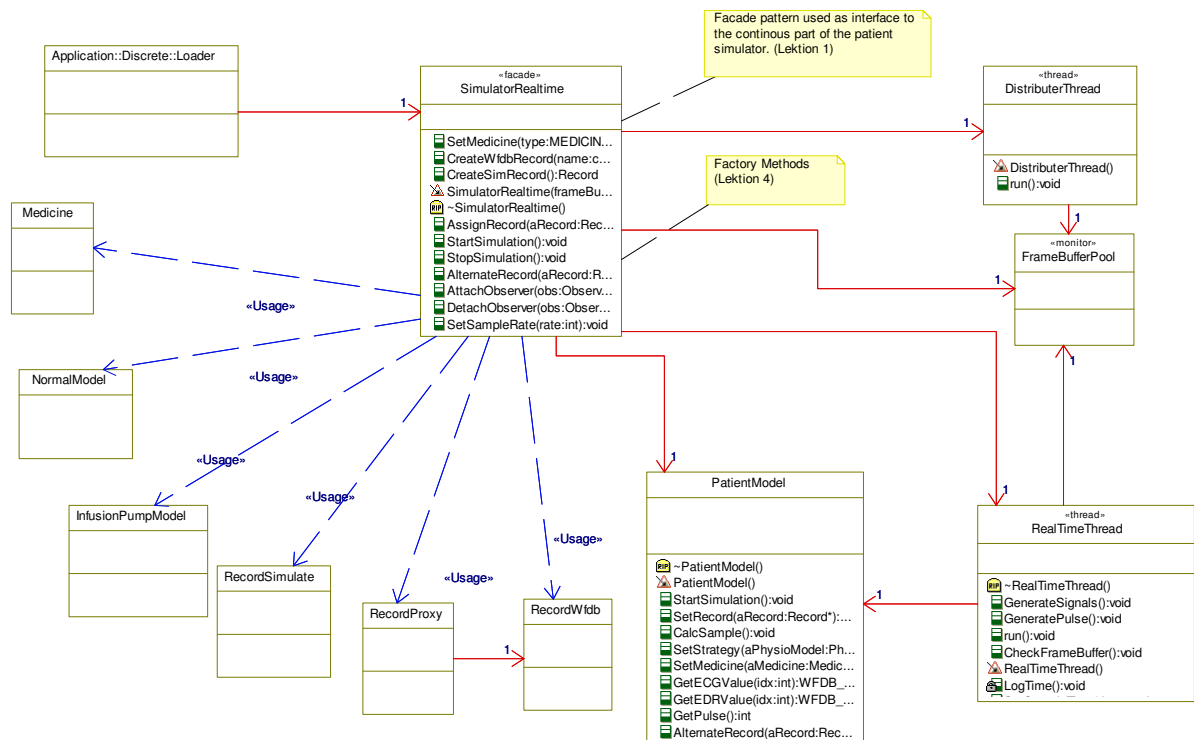


Figure 8 Façade Pattern used for interface to the Continuous Package

5.2.2 Discrete

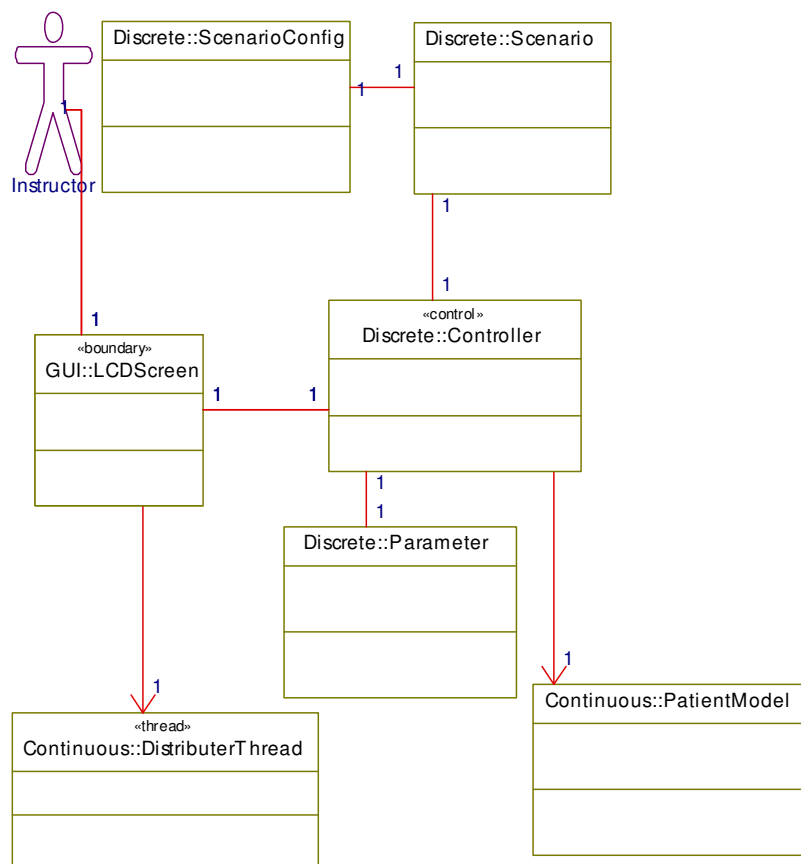


Figure 9 Application model for discrete package

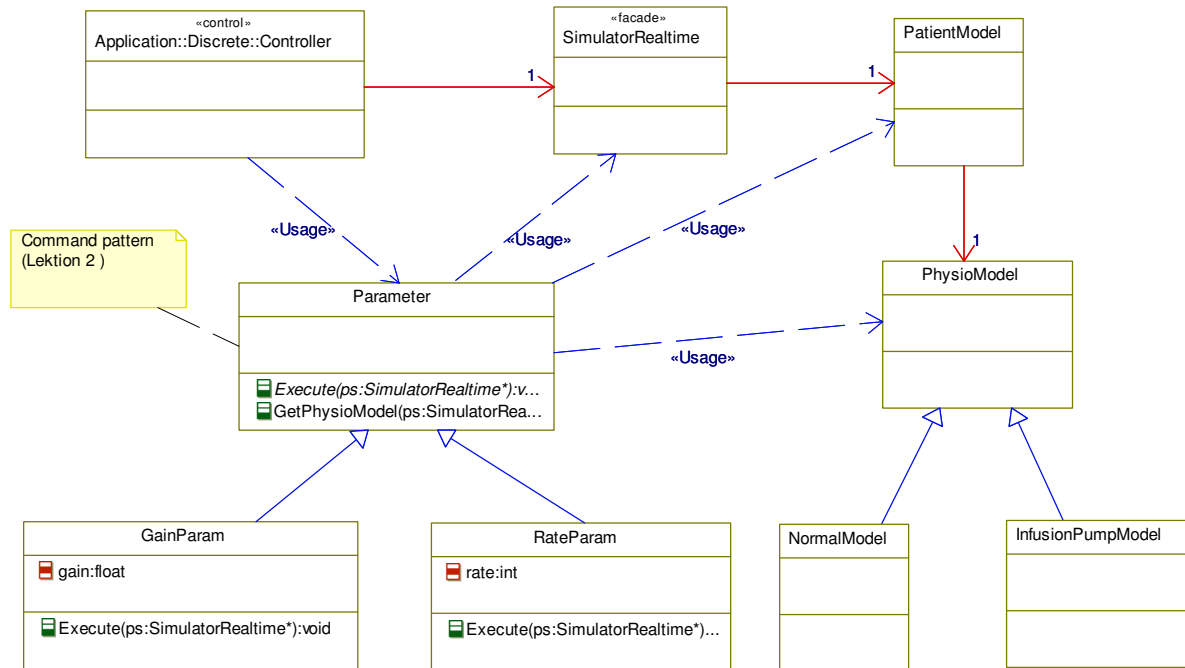


Figure 10 Command pattern used to set parameters for real time simulator

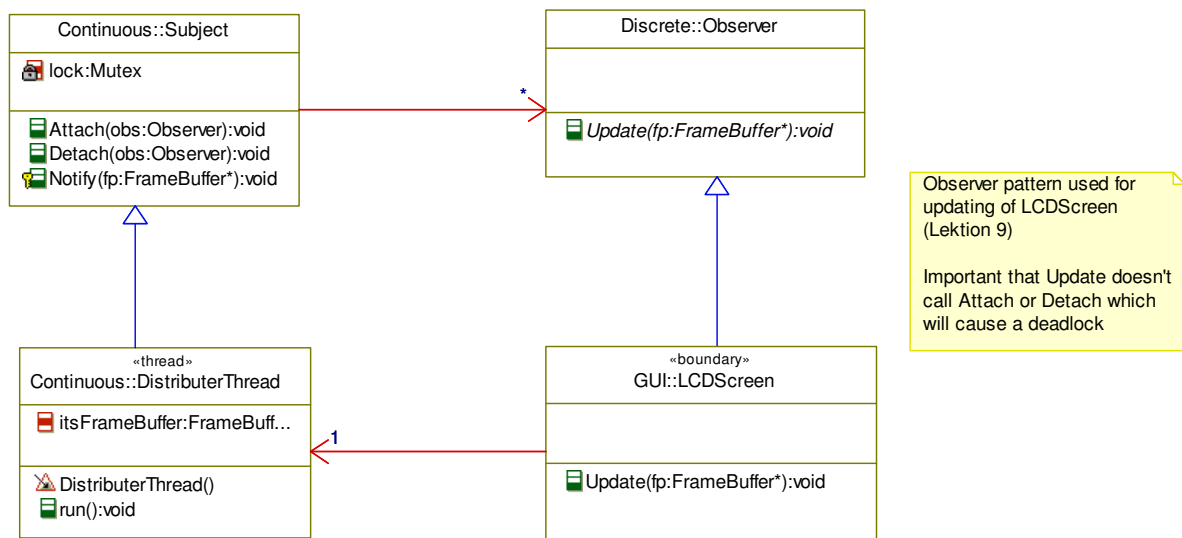


Figure 11 Observer pattern used to notify GUI with new framebuffer

5.2.3 Com

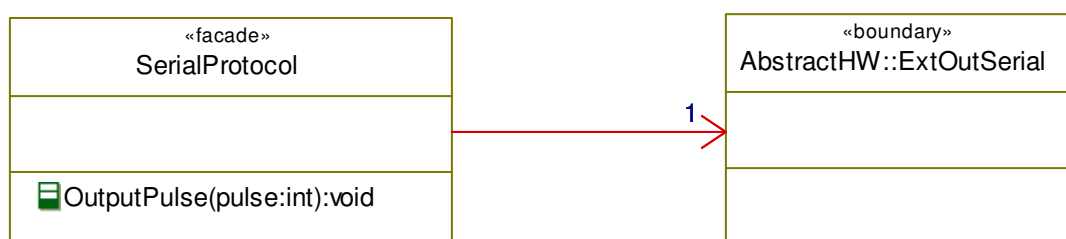


Figure 12 Communication package Serial protocol

5.2.4 AbstractHW

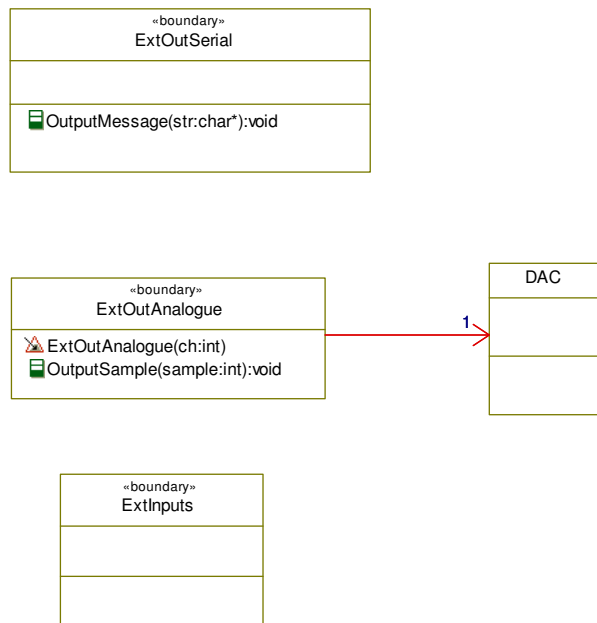


Figure 13 Hardware Abstraction

5.2.5 AbstractOS

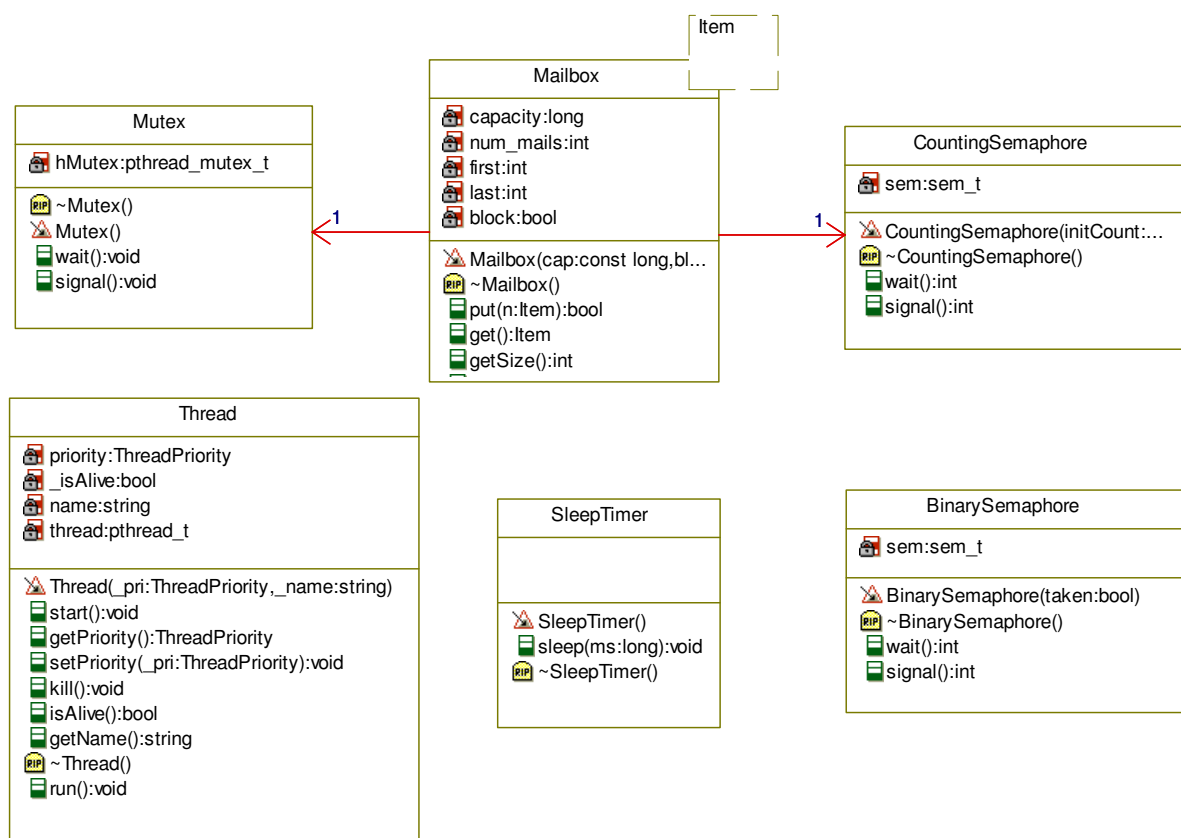


Figure 14 Abstract OS (Linux)

5.2.6 Application Helper Classes

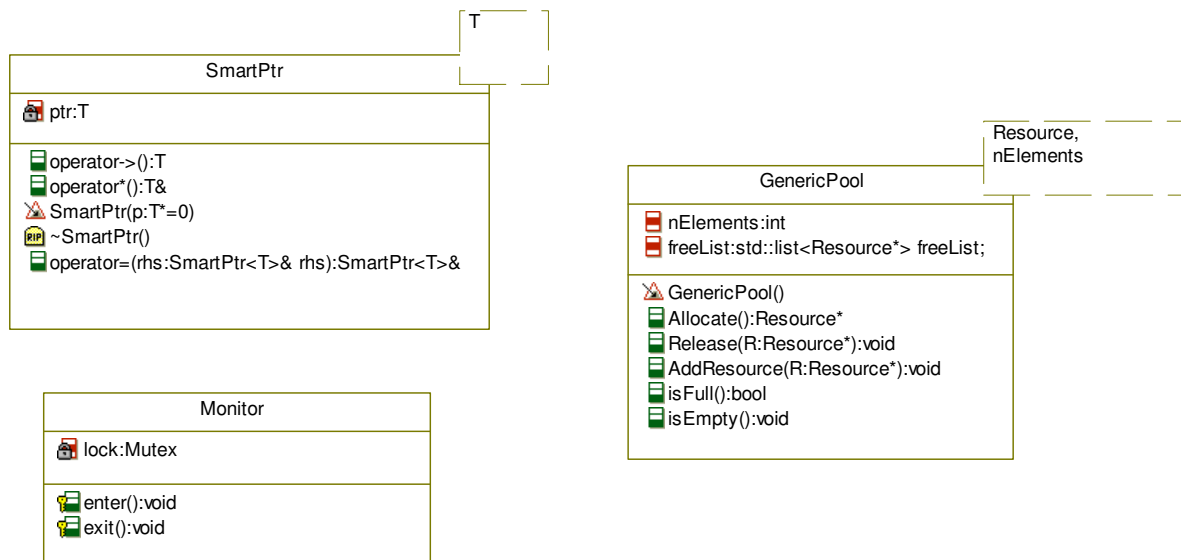


Figure 15 Application Helper classes

5.3 Use case realizations

5.3.1 Use case #1: Execute and Control Simulation scenarios

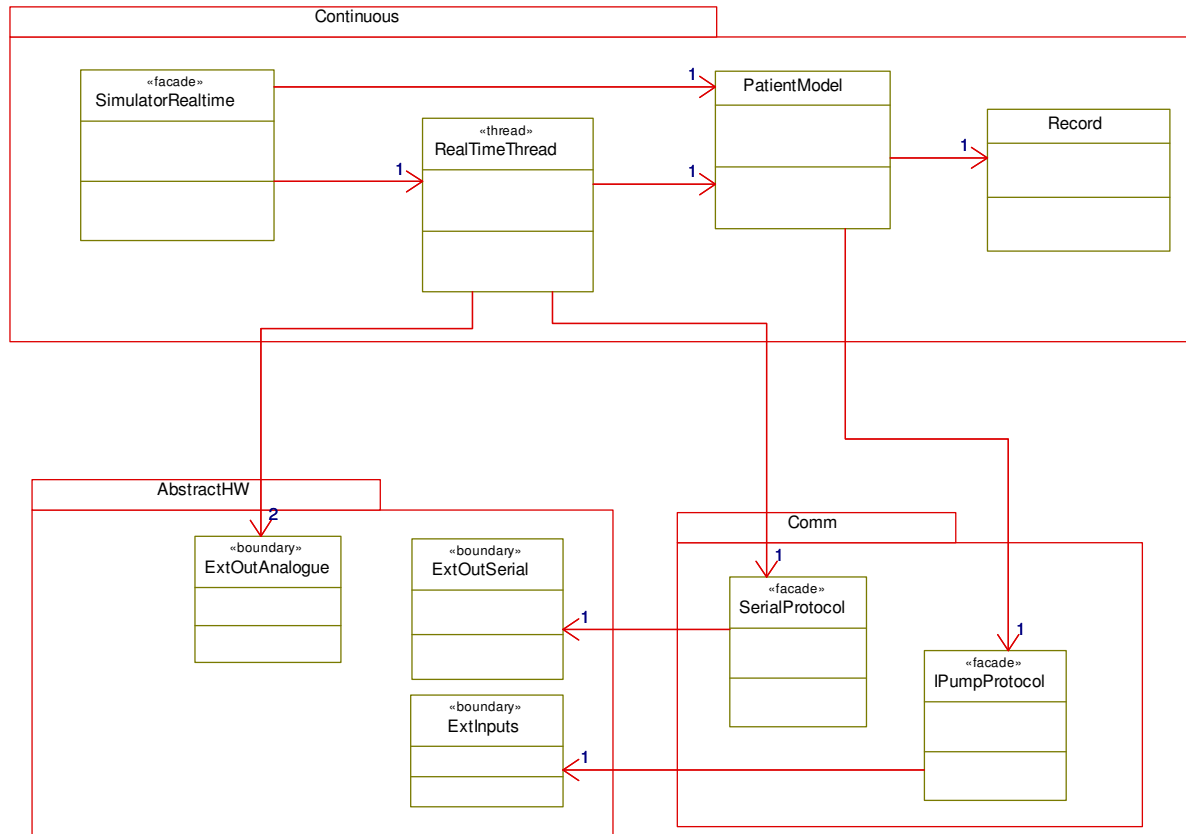


Figure 16 Logical view for use case #1 Execute and Control Simulation

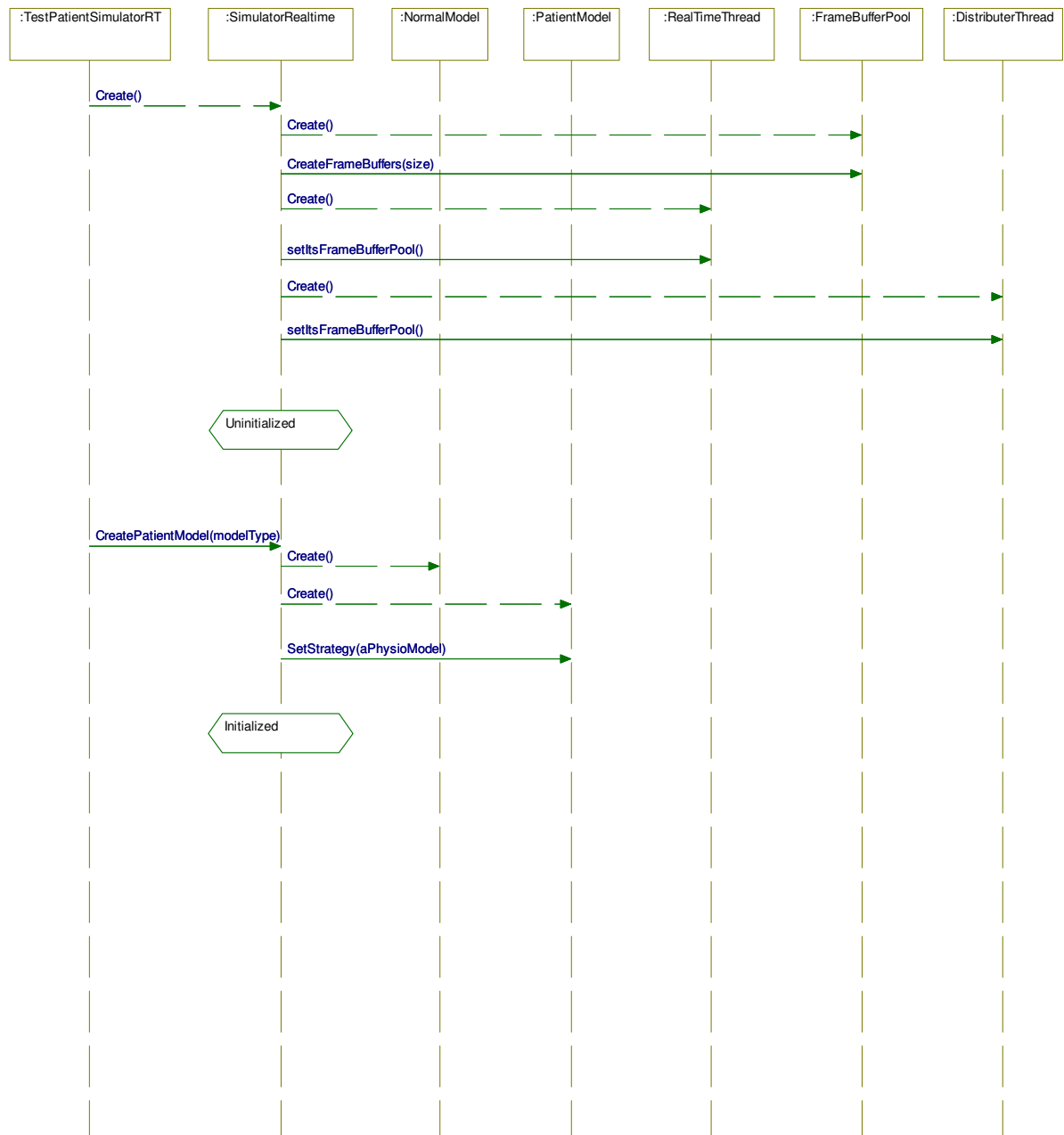
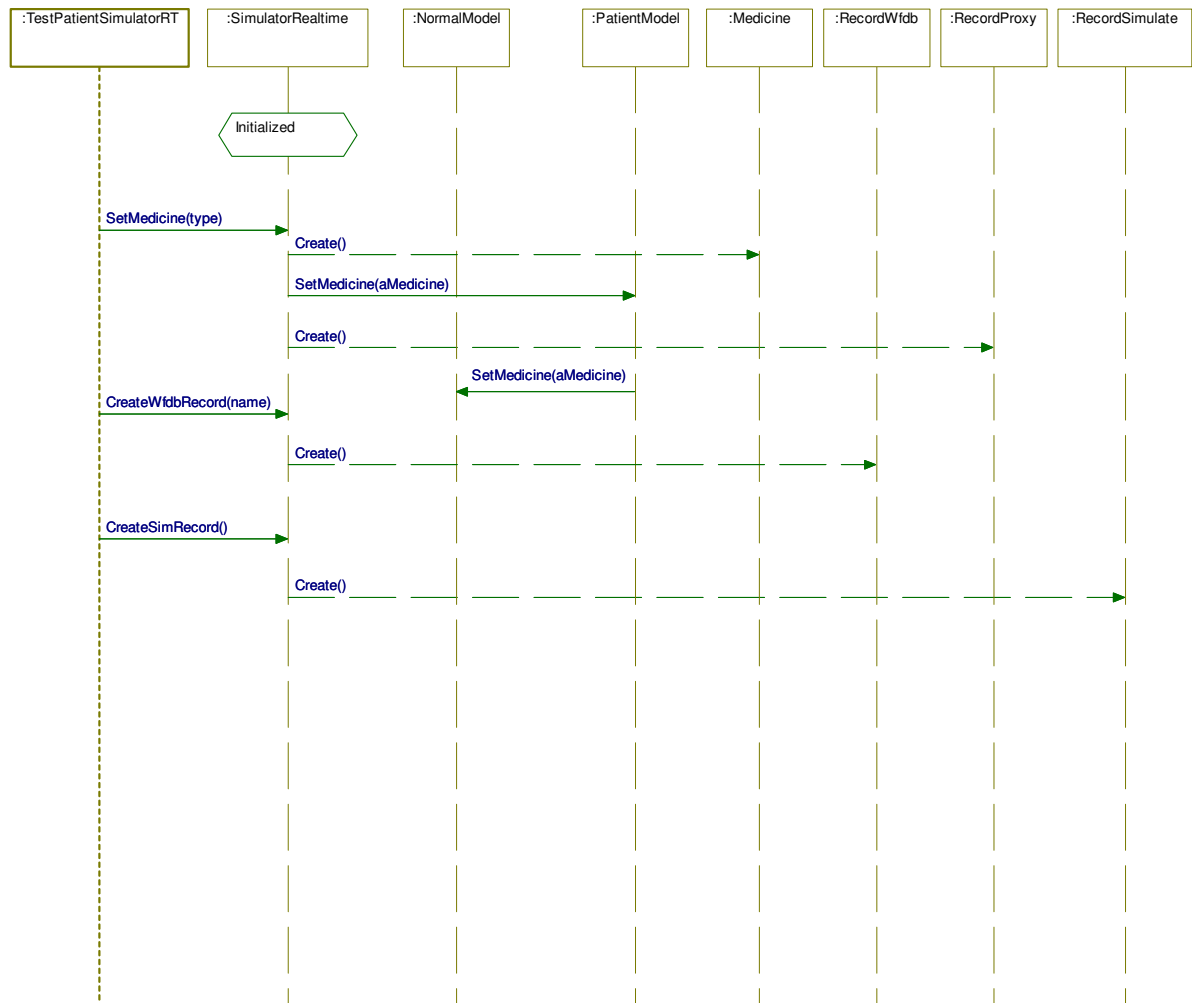


Figure 17 Sequence Diagram: Create Realtime Simulator

**Figure 18 Sequence Diagram: Create Simulation Record and Set Medicine**

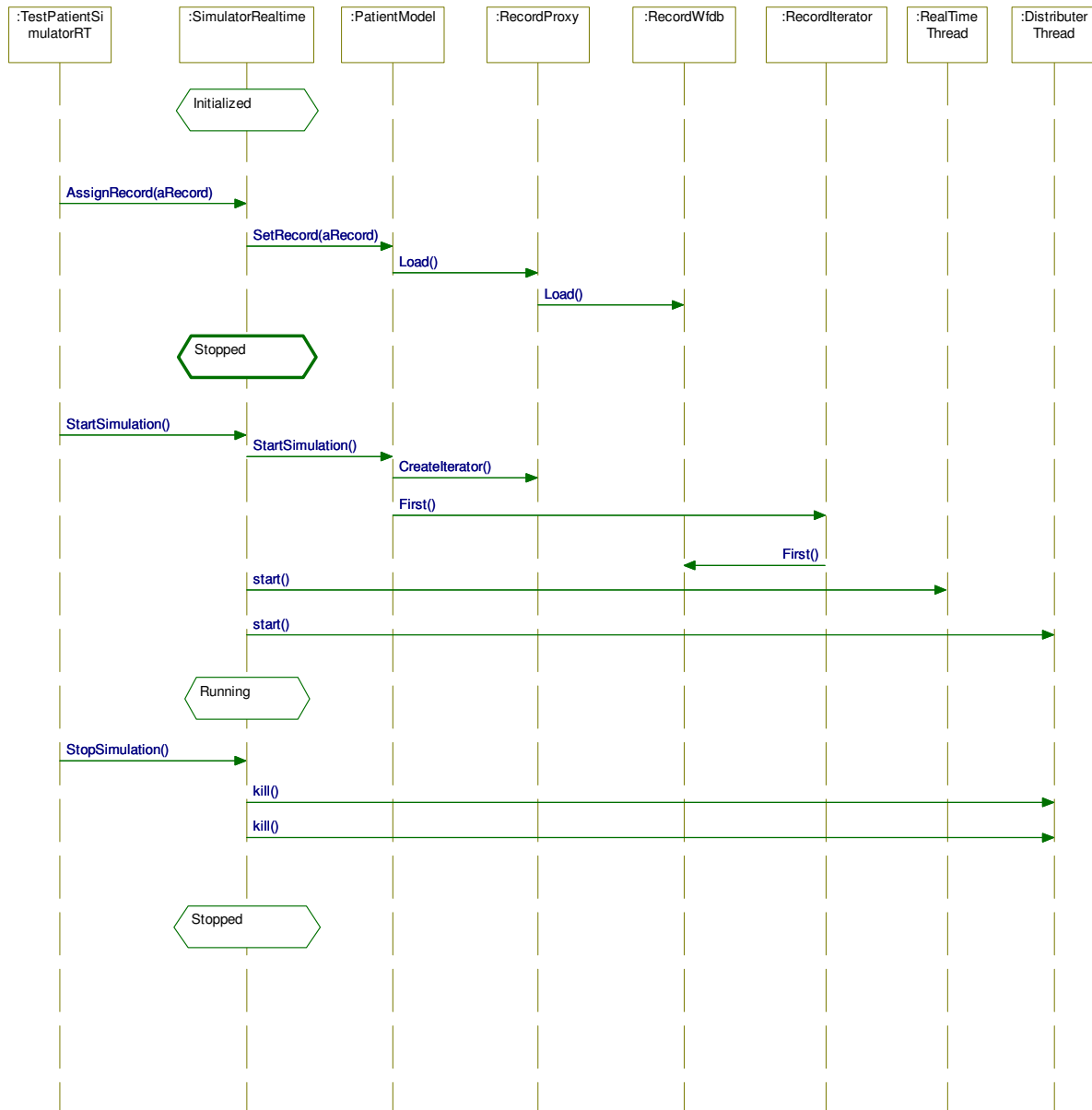


Figure 19 Sequence Diagram: Start Realtime Simulation

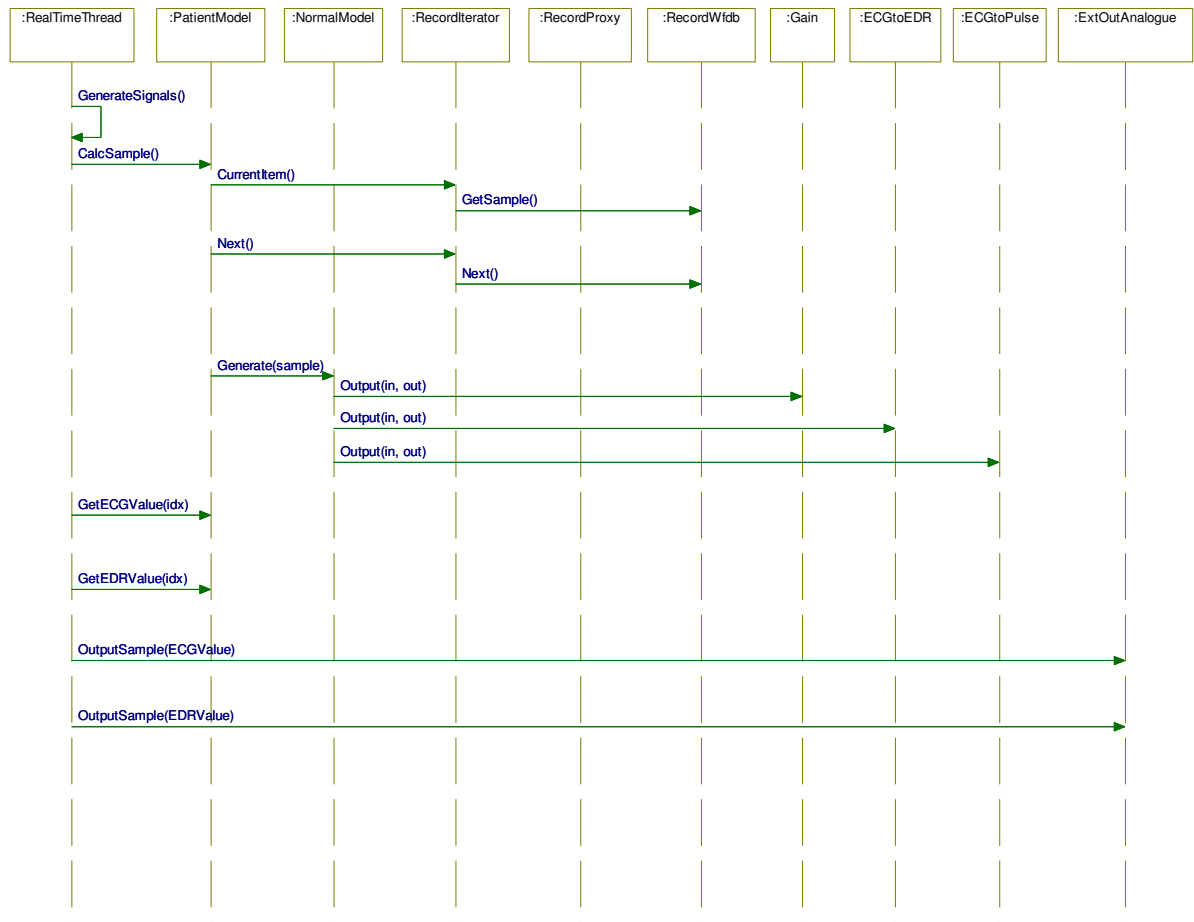


Figure 20 Sequence Diagram: RealTime thread GenerateSignals

5.3.1.1 Test setup for #UC 1 continuous package SimulatorRealtime

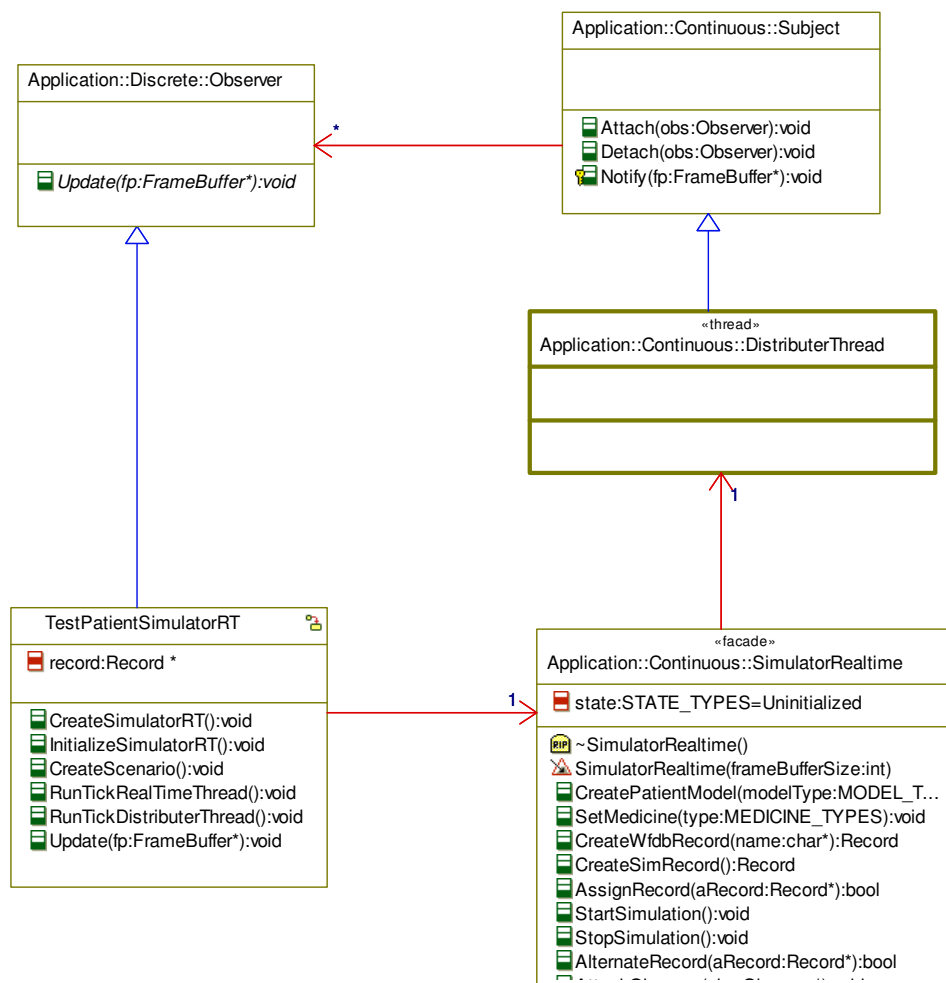


Figure 21 Test setup for UC#1 – test setup in Rhapsody only

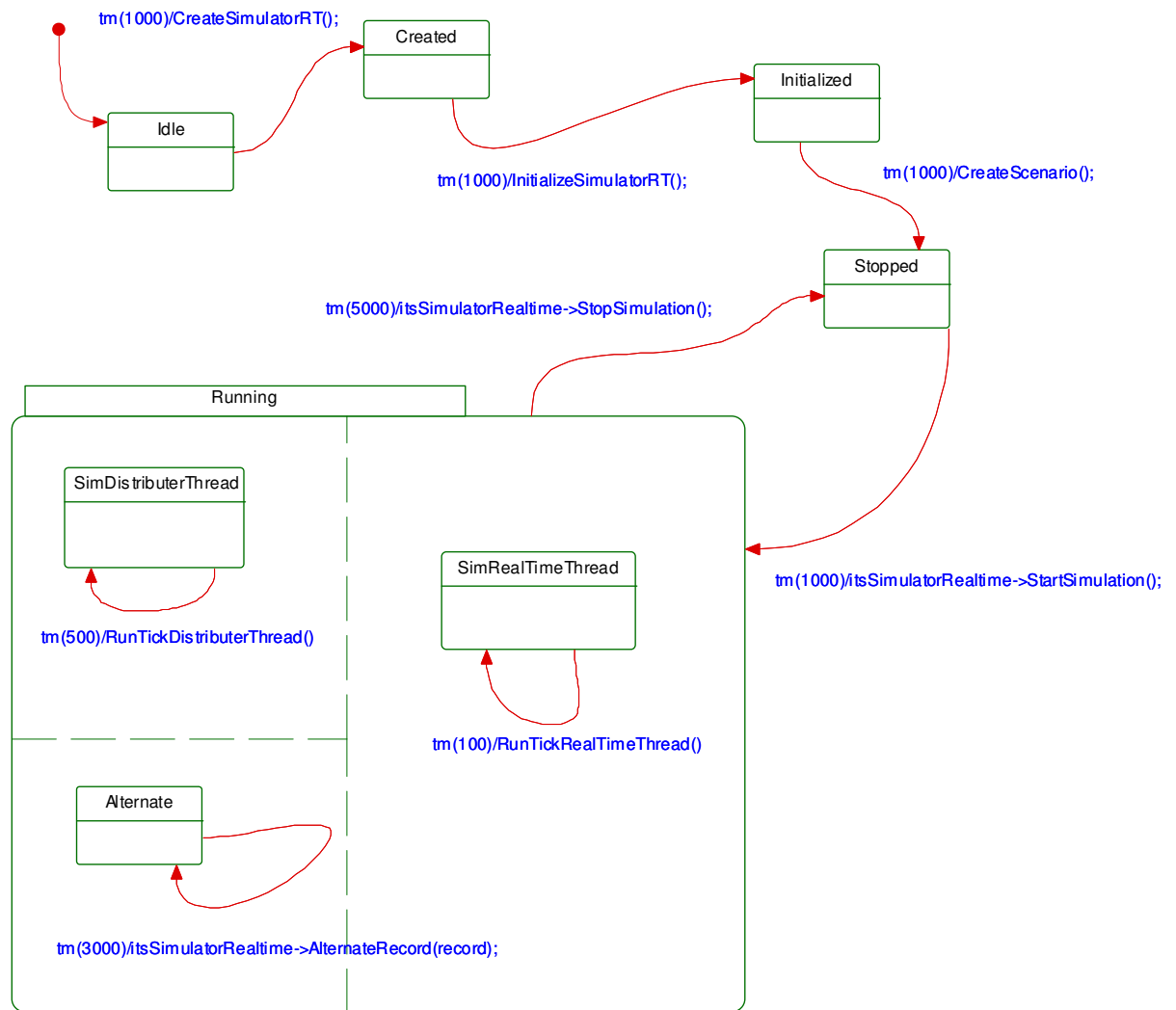


Figure 22 State machine test of UC#1 including simulation of Threads in Rhapsody

5.3.2 Use case 2. realization

6. PROCESS/TASK VIEW

6.1 Process/task overview

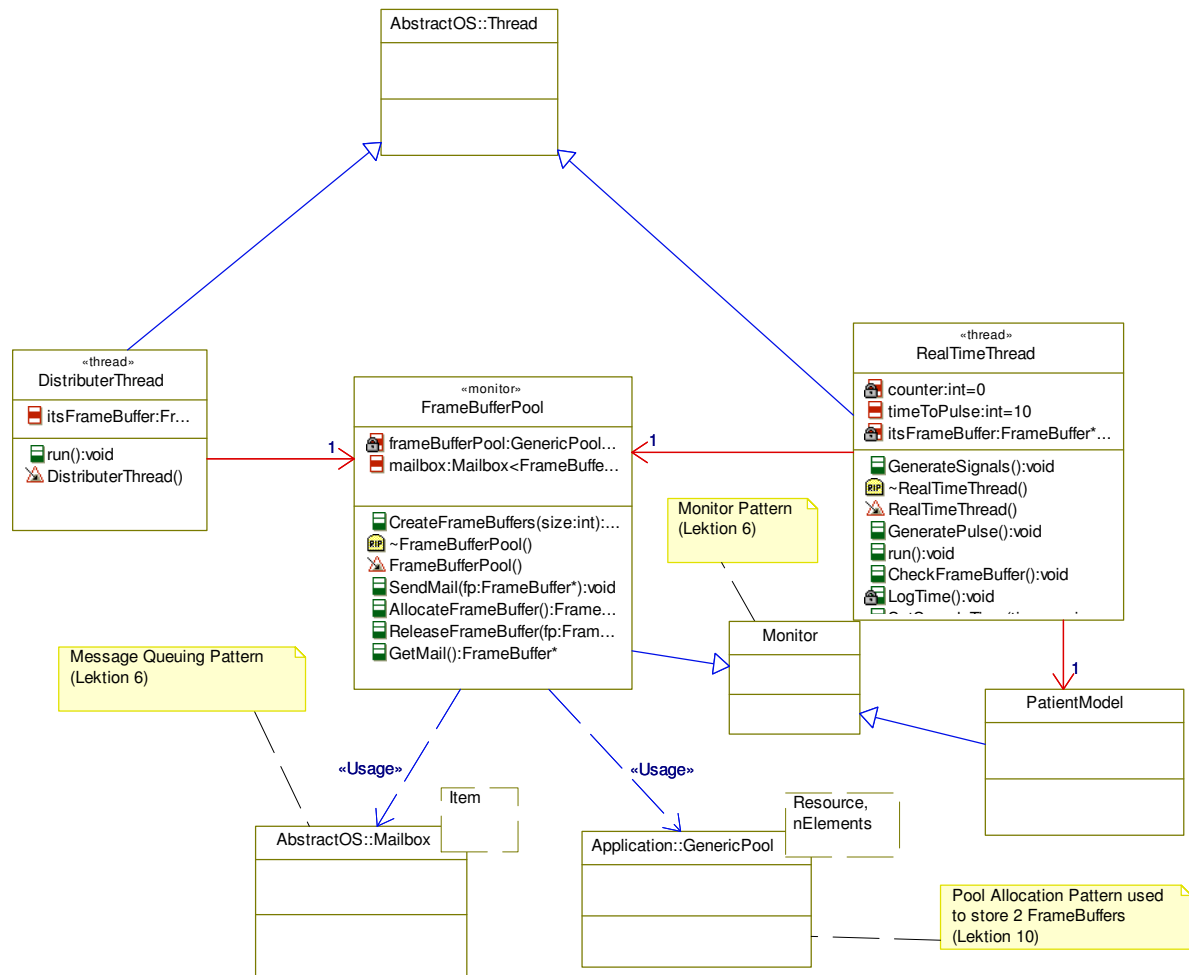


Figure 23 Process view for Distributer and RealTime threads and mechanism for synchronization

6.2 Process/task implementation

6.3 Process/task communication and synchronization

6.4 Process group 1.

6.4.1 Process communication in group 1


6.4.2 Process 1. description

6.4.3 Process 2. description

6.5 Process group 2.

7. DEPLOYMENT VIEW

7.1 System configurations overview

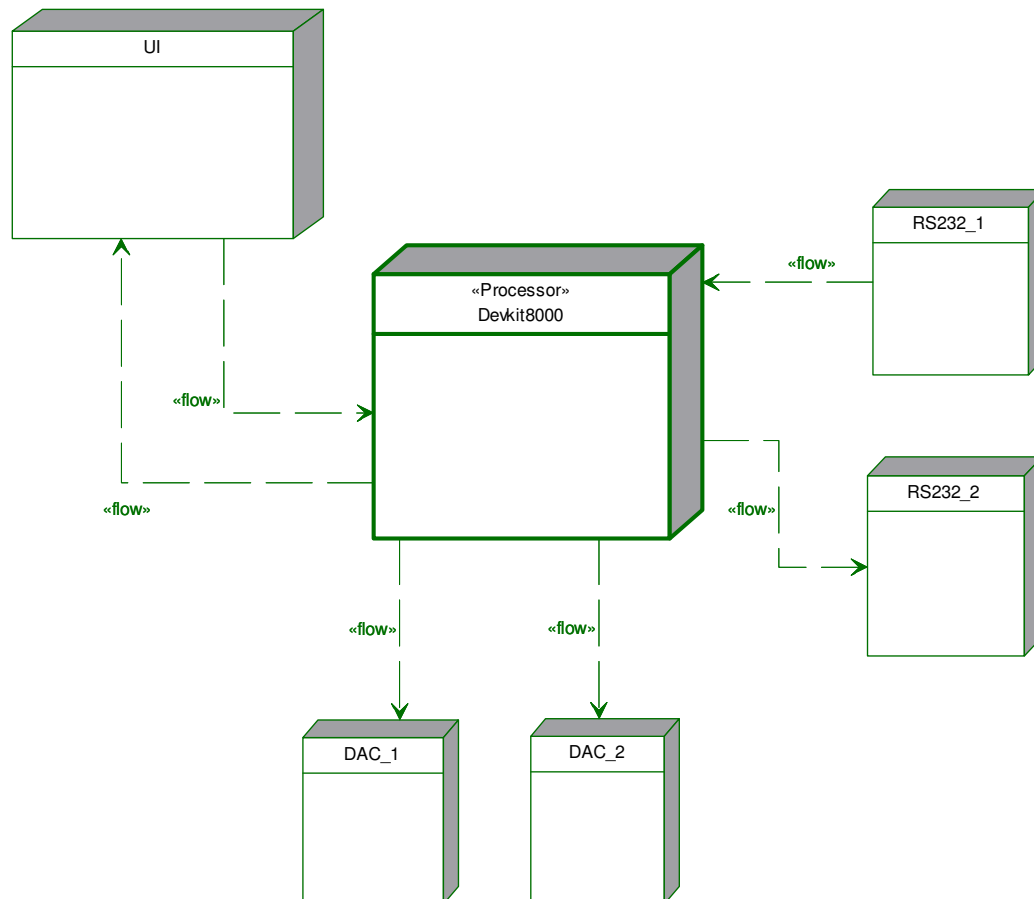

7.2 System configurations

7.2.1 Configuration 1.

7.2.2 Configuration 2.

7.3 Node descriptions

7.3.1 Node 1. description

7.3.2 Node 2. description

8. IMPLEMENTATION VIEW

8.1 Overview

8.2 Component descriptions

8.2.1 Component 1

8.2.2 Component 2

9. DATA VIEW

9.1 Data model

9.2 Implementation of persistence

10. GENEREL DESIGN DECISIONS

10.1 Architectural goals and constraints

The architectural goals are to make the system modifiable and provide high performance. Since the system is a real time system, there are constraints on how modular system can be, since modular and performance work in the opposite direction. It should be easy to extend the system new types of medicine, different kind of inputs and output. The system should be encapsulated so it's easy to port the system to different platform.

10.2 Architectural patterns

Observer Pattern:

Observer pattern help out with mass distribution of information in a system. If a number of processes are interested in the same information, observer pattern is the pattern to choose. In Sapien the observer pattern is used to distribute information generated by our simulation to the user interface. This a very common way to use observer pattern since it also notify the user interface when new data has arrived so it can be refreshed.

Five Layered Architecture:

Layered pattern organizes domains into a hierarchical organisation based their level of abstraction. The advance of layered architectural is that it make it easier to find relevant code, and make it easier to replace whole layers, for instance if you want to port the system to another device. There used open layered architectural so it's allowed to call more than one layer down.

10.3 General user interface design rules

10.4 Exception and error handling

10.5 Implementation languages and tools

Implementation languages:

- C++

Tools

- Eclipse
- TrollTech QT Creator

10.6 Implementation libraries

QT-Everywhere-4.6.1

wfdb

11. SIZE AND PERFORMANCE

12. QUALITY

13. COMPILATION AND LINKING

This project has been developed using a number of different tools.

IBM Rhapsody has been used for UML modelling, simulation and testing. It has been used for automatic source code generation and compilation on windows using Cygwin and generating C++ source code classes for Linux.

Eclipse has been used on Ubuntu Linux in VMware for implementation and testing of the WFDB interface and hardware interface to the DAC and serial ports.

Qt has been used on Ubuntu Linux for the final application by integrating and compilation of the graphical user interface, UML model and interface implementation made in eclipse.

The ARM cross compiler has been used to compile eclipse and Qt projects for the target.

In the following chapters the above tools usages for compilation the software for the Sapien patient simulator to target (Devkit8000) is described.

13.1 Rhapsody modeling and testing

Rhapsody (version 7.5) has been used to create an UML model for the Sapien patient simulator that is used for test of the simulator model including generated C++ source code before compilation on Linux and Target. We have created test scenarios by use of

the Rhapsody state diagrams and used them for testing the model by setting breakpoints in the states and made inspection of the state variables of the model by using the high level abstractions of that is possible with Rhapsody. This approach has reduced the development time in removing a lot of the manual C++ debugging and coding of the model. In this chapter we will describe how Rhasody has been configured for testing and code generation for Linux and Target.

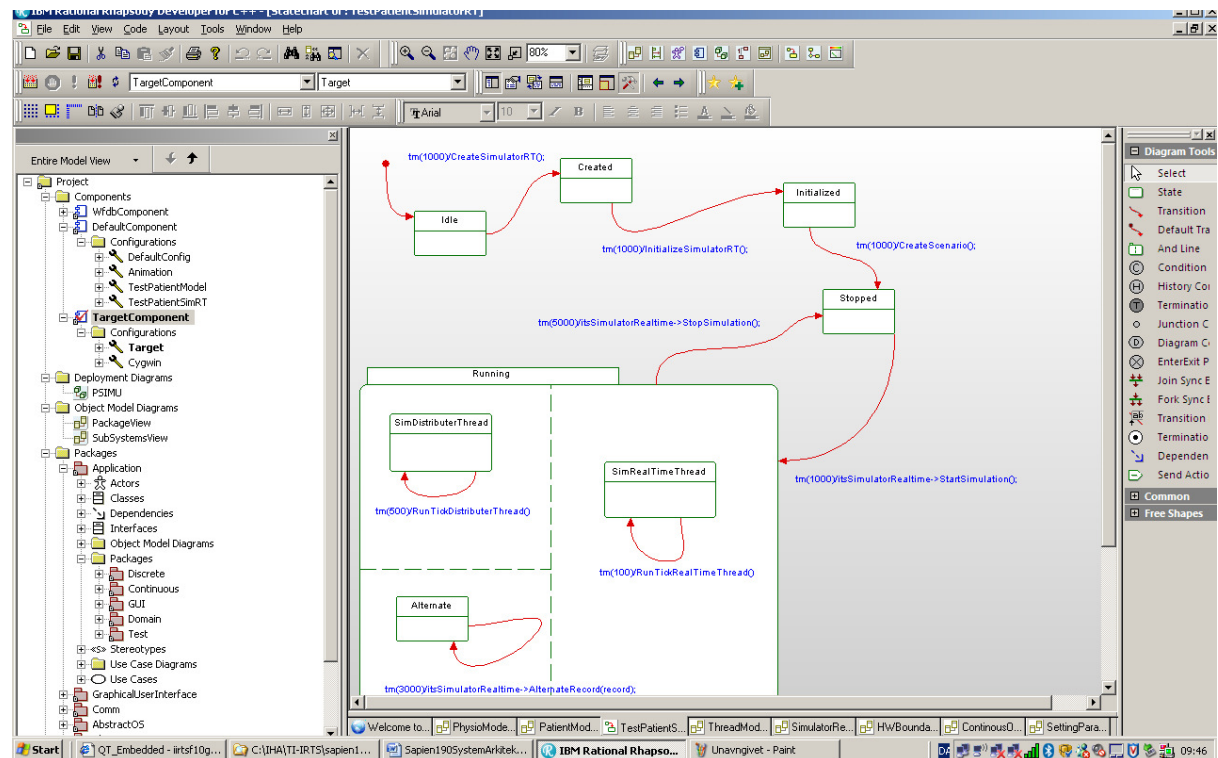


Figure 24 Rhapsody UML model for Sapien 190 used for simulation and test

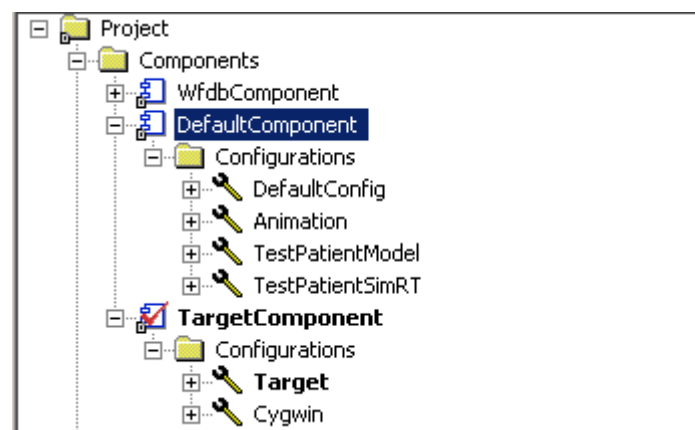


Figure 25 Rhapsody components for testing and source code generation

The figure above illustrates the components we have created for the Rhapsody project. The DefaultComponent and TargetComponent specifies the **scope** for classes for where source code is generated. Classes used for test of the model are not in scope of the TargetComponent since this component contains classes that are part of the final product. Simulation and test on windows is done by using Cygwin this is done by selecting the DefaultConfig settings to Cygwin and specifying the wfdb headerfiles that we have used

for the project:

DefaultConfig settings for windows simulation of DefaultComponent:

Standard Headers: wfdb/wfdb.h, math.h, wfdb/ecgcodes.h

Instrumentation Mode: None or Animation

Settings -> Environment: Cygwin

Rhapsody automatically generates the Makefile and compiles the generated source files in this configuration by using cygwin. We have made a configuration for each test scenario we have done in Rhapsody. (TestPatientModel, TestPatientSimRT)

TargetComponent is settings for generating source code for Linux and target:

Directory: C:\Ubuntu_share\sapient190\source\Sandbox\sapient_v1\cpy

Settings -> Environment: Environment

Here we have configured Rhapsody to generate the C++ source code directly to a sub directory of the Qt project that is compiled on Linux. We have shared a directory between the Windows and Linux platforms. The generated source classes is included in the Qt make project by using qmake.

13.2 Linux host Compilation-software

On the host Linux computer Qt and Eclipse must be installed. Qt has been used for the final product and Eclipse to test and develop parts of the Sapient software like hardware access and reading of WFDB records.

We have used the open source Qt Creator version 1.3.1 and the Qt SDK 2010.02 to be downloaded and installed from:

<http://qt.nokia.com/downloads>

We have used the open source Eclipse Galileo release to be downloaded and installed from:

<http://www.eclipse.org/cdt/downloads.php>

13.3 Linux Cross Compilation and linking process

The guide “Getting started with Qt” [9] describes how to install Qt on the Linux host which must be performed prior to the below steps. The following chapters describe how to cross compile the needed libraries for Qt and the WFDB library for the target Devkit8000 platform. For touch screen support see “Getting started with Qt”.

13.3.1 Qt Cross Compilation with qt-everywhere

This chapter describes how to download the embedded version of Qt that does not use the Linux X11 graphic library. It is therefore a suitable compact cross platform framework for the simulator platform where Devkit8000 is used for the first version of the product.

The following process describes how to configure Qt everywhere for compiling and linking a given program. Include a specification of special compiler and link switches.

This section is about how to modify qt-everywhere for cross compilation to target.

```
#
# qmake configuration for building for ARMv7 devices with arm-none-linux-gnueabi-g++
#
```

```
include(../../common/g++.conf)
include(../../common/linux.conf)
include(../../common/qws.conf)
```

modifications to g++.conf

```
QMAKE_CC = /opt/CodeSourcery/Sourcery_G++_Lite_2007q3/bin/arm-none-linux-gnueabi-gcc
QMAKE_CXX = /opt/CodeSourcery/Sourcery_G++_Lite_2007q3/bin/arm-none-linux-gnueabi-g++
QMAKE_LINK = /opt/CodeSourcery/Sourcery_G++_Lite_2007q3/bin/arm-none-linux-gnueabi-g++
QMAKE_LINK_SHLIB = /opt/CodeSourcery/Sourcery_G++_Lite_2007q3/bin/arm-none-linux-gnueabi-g++
QMAKE_CFLAGS += -O3 -march=armv7-a -mtune=cortex-a8 -mfpu=neon -mfloat-abi=softfp
QMAKE_CXXFLAGS += -O3 -march=armv7-a -mtune=cortex-a8 -mfpu=neon -mfloat-abi=softfp
#-mfpu=vfp
```

modifications to linux.conf

```
QMAKE_AR = /opt/CodeSourcery/Sourcery_G++_Lite_2007q3/bin/arm-none-linux-gnueabi-ar cqs
QMAKE_OBJCOPY = /opt/CodeSourcery/Sourcery_G++_Lite_2007q3/bin/arm-none-linux-gnueabi-objcopy
QMAKE_STRIP = /opt/CodeSourcery/Sourcery_G++_Lite_2007q3/bin/arm-none-linux-gnueabi-strip

QMAKE_INCDIR += /home/stud/tslib_arm/include
QMAKE_LIBDIR += /home/stud/tslib_arm/lib
```

```
load(qt_config)
```

Configure the libraries to use the linux-armv7-g++ configuration:

```
$ cd /home/stud/qt-everywhere-opensource-src-4.6.2
$ ./configure -embedded arm -xplatform qws/linux-armv7-g++ -qt-kbd-linuxinput -qt-mouse-
tslib-opensource -verbose -R /home/stud/tslib_arm/lib/
```

Make the libraries (Grab ...something, this takes several hours):

```
$ make
```

and install them:

```
$ make install
```

13.3.2 WFDB Installation and Cross Compilation

This chapter describes how to install the WFDB library from physionet.org and how to modify it for cross compilation to the ARM target running Linux. The WFDB library is an open source project used to read and manipulate patient records from the PhysioBank database. The host Linux platform must have installed the “arm-none-linux-gnueabi-gcc” prior to this installation.

First download the WFDB library from here:

<http://www.physionet.org/physiotools/wfdb.shtml#downloading>

Steps to install WFDB on the Linux host:

1. Install XView for Linux
\$ sudo apt-get install xviewg xviewg-dev
2. Unpack the downloaded wfdb.tar.gz
/home/stud\$ tar -xvf wfdb.tar.gz
3. Run configure in the wfdb directory
/home/stud/wfdb-10.5.1\$./configure
4. Build and make the wfdb library
/home/stud/wfdb-10.5.1\$ make
5. Install the wfdb library on Linux host
/home/stud/wfdb-10.5.1\$ sudo make install

Steps for cross compile of the WFDB library to target:

This part describes how create the libwfdb.so.10.5 library that must be copied to the target platform. Prior to this step the WFDB library must be installed on the Linux host see description above.

1. The first step is to copy the source library files to a new directory. Copy the lib directory from the wfdb-10.5.1 installation to a new directory named lib-arm.

```
/home/stud/wfdb-10.5.1/lib
Copy to
/home/stud/wfdb-10.5.1/lib-arm
```

2. Edit the Makefile manual in the new lib-arm directory, change the following lines to:

```
SRCDIR = "/home/stud/wfdb-10.5.1"
WFDBROOT = /home/stud/wfdb-arm
CC = arm-none-linux-gnueabi-gcc
BUILDLIB = arm-none-linux-gnueabi-gcc $(MFLAGS) -shared -Wl,-
soname,$(WFDBLIB_SONAME) $(LL) \
-o $(WFDBLIB)
```

These modification will now used the arm cross compiler (arm-none-linux-gnueabi-gcc) and install the WFDB header and library files in /home/stud/wfdb-arm

3. Make, Compile, link and install the wfdb library for arm target

```
/home/stud/wfdb-10.5.1/lib-arm $ make
/home/stud/wfdb-10.5.1/lib-arm $ make install
```

Ignore errors when “make install” is performed this is due to some files it tries to install on the Linux host that is not needed on target.

Steps for using the WFDB library in eclipse:

1. Change the eclipse application project to include wfdb:
/home/stud/wfdb-arm/include
2. Change the eclipse project to include the wfdb library and path:
-l wfdb
/home/stud/wfdb-arm/lib

Steps for using the WFDB library in Qt:

This step describes how to use the qt-everywhere version of qmake to generate the makefile for cross compilation to target. The following parameters to qmake must be added as described below.

These parameters adds information about the paths for find the wfdb header files and wfdb library. The below line also specifies to include the touch screen library (-lts) for the Devkit8000 platform.

```
/home/stud/qt-everywhere-opensource-src-4.6.1/bin/qmake LIBS+="-L/home/stud/wfdb-arm/lib -lts -lwfdb" DEFINES+=_LINUX DEFINES+=_USE_HW_DAC
INCPATH+=/home/stud/wfdb-arm/include
```

If the _USE_HW_DAC is not specified a version will be compiled that just prints out record samples on the standard output. (Terminal output)

14. INSTALLATION AND EXECUTING

This section describes how to install the sapien190 software on the DevKit8000 target running Linux. We will describe how to install the needed libraries and the application program on target.

14.1 Installation

The sapien190 executable can be copied to the target by copy of the cross compiled binary executable program to the SD card or using Ethernet over an USB connection to the target.

Ethernet over USB to target

Configure Ethernet over USB after the USB wire is connected between host and target powered on with the target IP address of 10.9.8.1:

First configure the IP address of the Linux host:

```
$ ifconfig usb0 10.9.8.2/24 up
```

Use secure copy of binary cross compiled application to target in the default directory /home/root:

```
$ scp sapien190 root@10.9.8.1:
```

Copy of files using SD card:

Copy the sapien190 binary file to the SD card and insert it on the target. The SD card will automatically be mounted on the target at:

```
/media/mmcblk0p1/
```

14.2 Executing-hardware

To run sapien190 on the host development platform Qt and the WFDB library must be installed to start sapien190 enter the following command in Linux:

```
$/sapien190
```

To run Sapien190 on the target after installed:

```
$/sapien190 -qws
```

14.3 Executing-software

This chapter describes the needed libraries to be installed on target to run the Sapien190 patient simulator.

Installing Qt Everywhere on target by moving the cross compiled libraries and some fonts to the DevKit8000 target.

Fonts: At least one font must be present on target to show text. Fonts can be found in `$HOME/qteverywhere-opensource-src-4.6.1/lib/fonts` on the host and must be copied to `/usr/local/Trolltech/QtEmbedded-4.6.1-arm/lib/fonts/` on target (create the directories as necessary).

(At least) the below libraries must be copied from `$HOME/qt-everywhere-opensource-src-4.6.1/lib/` on host to `/usr/lib` on target:

- `libQtCore.so.4`
- `libQtGui.so.4`
- `libQtNetwork.so.4`

The wfdb library must be copied from host `/home/stud/wfdb-arm/lib/` on host to `/usr/lib` on target:

- `libwfdb.so.10`

The standard C++ library must be copied from host `/home/stud/wfdb-arm/lib/` to `/usr/lib` on target:

- `libstdc++.so.6`

Install and enable the touch screen on Devkit8000 by following the steps in “Getting started with Qt” [9] and install the touch screen library as described below:

Copy the touch screen libraries from `/home/stud/tslib_arm/lib/` to `/usr/lib` on target and enable the touch on the DevKit8000:

- `libts-0.0.so.0`
- `libts.so`

Enable touch screen on DevKit8000:

```
$ chmod a+rw /dev/input/event2
$ export QWS_MOUSE_PROTO=Tslib:/dev/input/touchscreen0
```

Add the environment variable `QWS_MOUSE_PROTO` in the startup profile listed below.

```
/etc/profile
```

14.4 Execution-control (start, stop and restart)

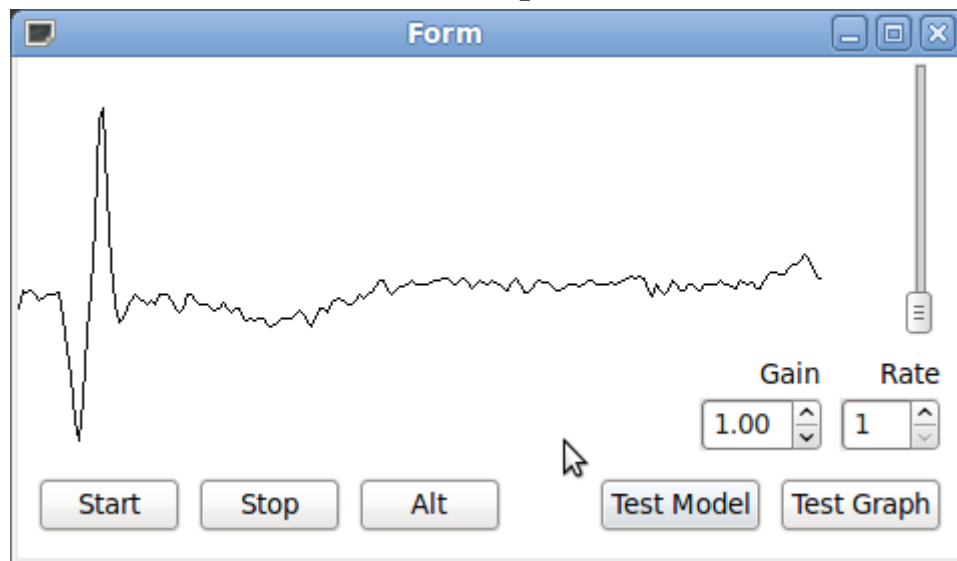


Figure 26 Sapien190 main dialogue

The first version of the Sapien190 patient monitor is only able to start and stop playing a specific patient ECG record e0104 alternating with patient ECG record e0103. When test model or test graph is selected the contents of patient record files is printed on the display independent of the DAC signal outputs. The gain and rate of playing the record file can be adjusted during simulation.

14.5 Error messages

15. APPENDICES